



# ServiceNet Medical X-ray

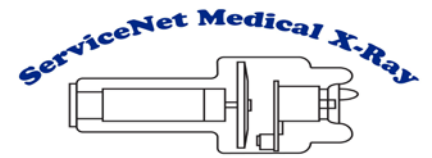

**ViVIX-S** is a 17x17 flat panel using a 9 Mega-pixel active digital radiographic area and 14-bit signal digitization. The **ViVIX-S** system ensure that with a large field coverage area and software designed for Veterinary radiographic applications with the use of its unique image processing system and proprietary flat panel detector technology you get the best possible image the first time not the second.

**ViVIX-S** can be used without any connections to the X-ray generator with the patented AED (Automatic Exposure Detection) circuit included on all Vieworks flat panel detectors. **ViVIX-S** works just like CR or film as it is triggered by the radiation and not by a generator interface connection.

User friendly Veterinary imaging software with DICOM 3.0 standard provides consistent image quality at a dramatically reduced dose and faster image information with optimized algorithms for each different study. Image can be simply acquired and transmitted to the DICOM server through the communication cable in seconds.

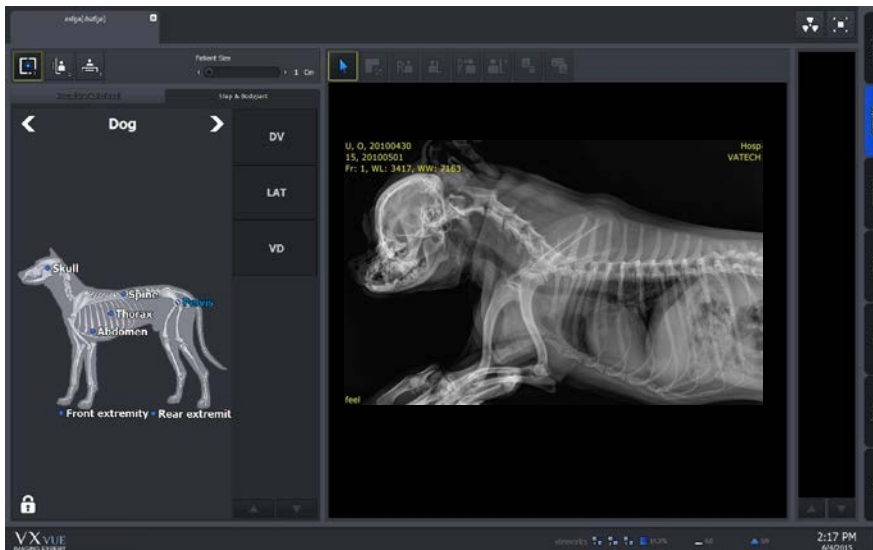

5262 Washington Rd. Albany, OH. 45769 (Telephone) 1-800-830-1035 (Fax) 1-866-887-9729

Email: [Rob@xrayservice.com](mailto:Rob@xrayservice.com)

Website: [www.servicenetmedical.com](http://www.servicenetmedical.com)